


## Precautions When Using CBLR-CCP20

Observe the following precautions when using CBLR-CCP20.

- To use CBLR-CCP20, set the image input board (frame grabber) to PoCL.
- Do not connect the Camera Link cable 2 m to the image input board that is set to PoCL.
- When using the Camera Link cable 2 m, set to OFF status the PoCL settings of the image input board.

\* Refer to “OLS4000 Safety guide” when setting the image input board.

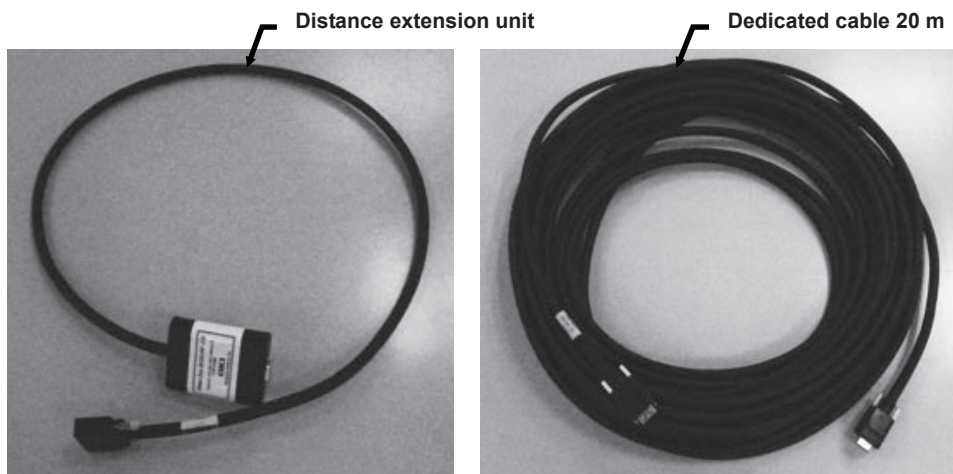

CBLR-CCP20



Camera Link cable 2m

(日本語は裏面をお読みください)



AX7986