

## U-P5REMC

Motorized Centerable Quintuple Revolving Nosepiece with Slider Slot for DIC



### SPECIFICATIONS

Global item code / Local code for USA		N2662800/U-P5REMC-1-5
Predecessor product information		None
Product Category		Nosepiece
Nosepiece	Method of rotating	Motorized rotation
	Quantity of the holes	5
	Centerable	Available (4 holes can be centered.)
	Slider slot	Available
	Encoded	5 holes
	Switching time	Within 0.65 s
	Power supply method	Supplied by Control box
Attached method	Illuminator or Arm	slide dovetail
	Objective lens	W 20.32 mm × 0.706 mm
	Slider	Slider size: 34 mm × 13.5 mm Length from the bottom surface of opening to mounting surface: 19 mm
	Cable connection	
	Control box(BX3M-CB)	Intermediate cable connection(BX3M-RMCBL: 0.2 m)
	Control box(BX3M-CBFM)	Intermediate cable connection(BX3M-RMECBL: 2.9 m)
ESD function		–
Weight		About 1.1 kg
Country of Origin		Japan

### SYSTEM ENVIRONMENT LIMITATIONS FOR OPERATION

Environment	Indoor use
Temperature T0 [°C]	5 – 40
Humidity RH0 [%]	Maximum 80% (up to 31°C) (no condensation) Usable humidity declines linearly as temperature rises above 31°C 34°C (70%) to 37°C (60%) to 40°C (50%)
Power supply and voltage stability	±10%
Pollution level (IEC60664)	2
Over voltage category (IEC60664)	II

# U-P5REMC

## RESTRICTIONS

The motorized nosepiece cannot be rotated manually.


When detaching the motorized nosepiece from the microscope frame in order to replace the objective, turn OFF the power of the microscope frame in advance.

Attach the objective protection cap to the nosepiece hole that is not in use.

## PACKING LIST

Item		Quantity
Main unit	U-P5REMC-B	1
Accessories	Centering wrenche	2
	Attention tag	1

## MAIN DIMENSION

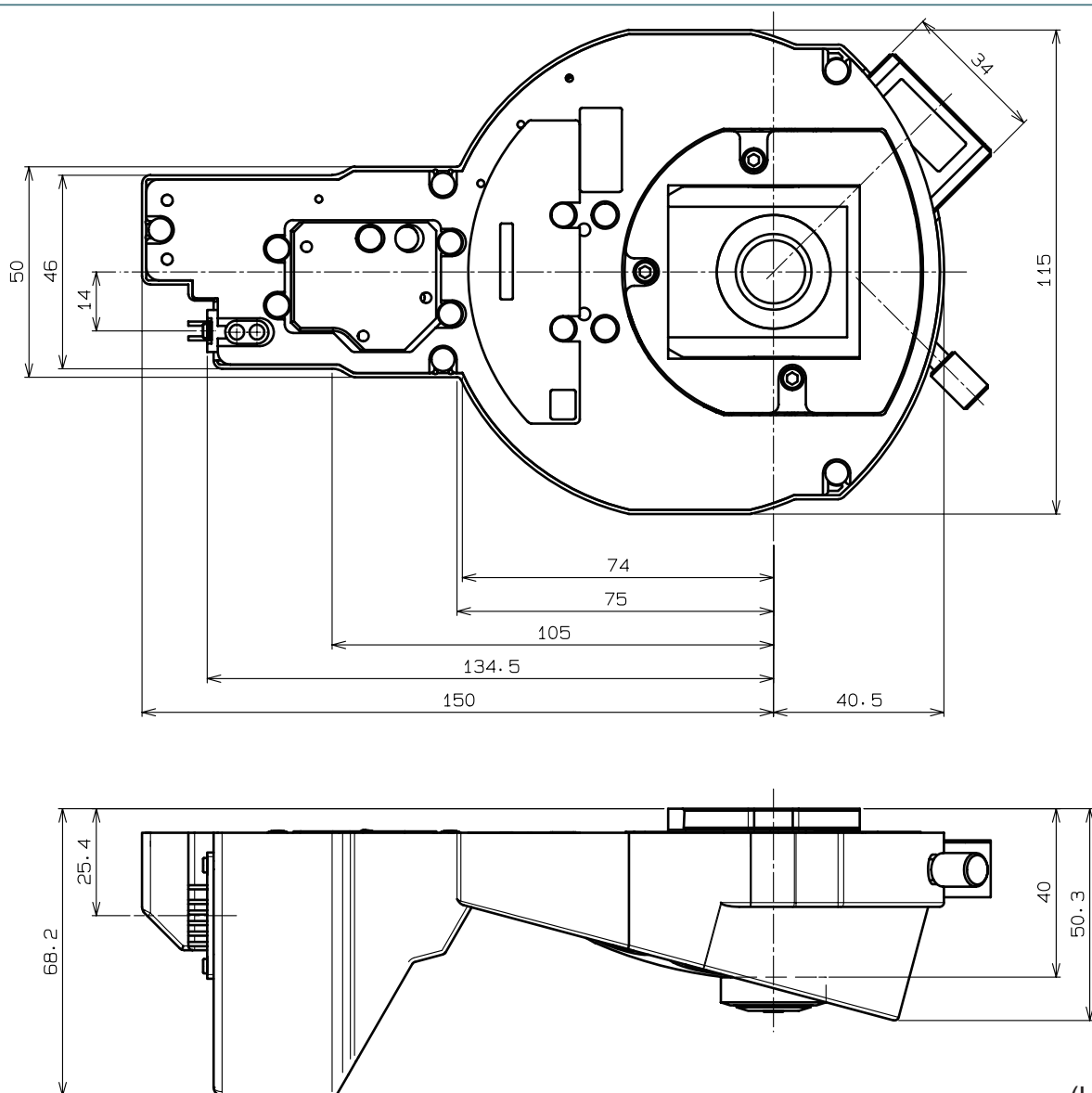

(Unit: mm)

· Specifications and appearances are subject to change without any notice or obligation on the part of the manufacturer.